

ELEKTRA™

— *comfort* —



INSTRUCTION MANUAL HEATING PAD



Read these instructions before use and keep them safe. If you pass the appliance on, pass on the instructions too. Remove all packaging, but keep it until you are satisfied.

Important Safeguards

Before connecting the appliance check that the voltage indicated on the label corresponds to the mains voltage in your country.

When using electrical appliances basic safety precautions should always be followed. Please read the following carefully before use:

1. After unpacking, make sure that there are no visible signs of damage.
2. This appliance must only be used as a heat pad as indicated in the instructions. Any other use is to be considered hazardous. Pharmalec Agencies C.C. and any of its associate companies cannot be held responsible for any damage or injury that may occur due to improper and incorrect use.
3. Do not use with wet hands or feet, or while your body is wet.
4. Do not use the Heat Pad if it is wet.
5. Do not pull the appliance or the power cord to remove the plug from the socket.
6. Before cleaning the appliance, disconnect it from the mains by removing the plug.
7. Should the Heat Pad or protection fuse fail to work properly, turn off the appliance and unplug it from the power outlet. Contact the Consumer Service Centre for examination and repair.
8. Always unplug the appliance immediately after use.
9. Check the cord regularly for wear or damage, particularly where the cord enters the appliance and the plug.


10. During use, ensure that the power cord is completely unwound.
11. Close supervision is necessary when this appliance is used by children and invalids.
12. Always allow the appliance to cool down before packing away.
13. Only use the Heat Pad with the controls provided, as it is part of the overheat protection system.
14. Do not fold the Heat Pad.
15. Avoid creasing the Heat Pad.
16. Do not iron.
17. Please follow cleaning instructions carefully so as to avoid injury or damage.
18. This appliance is not designed for medical or hospital use.
19. Do not insert needles or pins into the Heat Pad.

Household use only

Additional Safeguards

1. Do not try to repair the appliance yourself. Contact the approved Service Centre for examination and repair.
2. This product has been designed to operate safely when used according to the instruction guide.
3. If the electrical circuit is overloaded with other electrical appliances, your Heat Pad may not operate properly. It should then be operated on a separate circuit from other appliances.

Know your Elektra Comfort Heat Pad

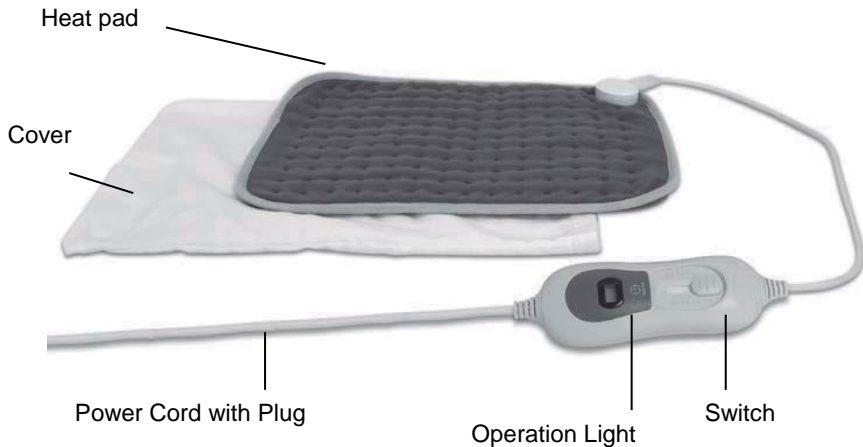

Turning on your Heat Pad and Temperature Selection



- To turn on the heat pad slide the button up from **0**.
- To select the heat you require, slide the button to either heat level **1**, **2** or **3** as desired.
- To turn off the heat pad slide the button down to **0**.
- When the heat pad is on the chosen setting will be displayed in the selection indicator screen.


N.B. Do not place the control switch underneath a pillow or a blanket when it is turned on. The control switch must never be covered.

Instructions for Cleaning

Heat pad cover:



Heat pad:

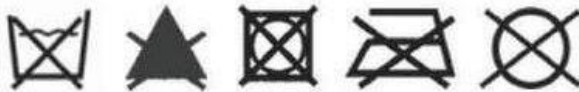

Should you need to clean the heat pad:

- Wipe it down with a damp cloth and a small amount of non-abrasive cleaning agent.
- Never get the heat pad wet.
- Do not wring it out
- Do not use it if it is wet.

Storage Instructions

- The best place to store your Heat Pad when not in use, is in the box you purchased it in.
- Should you wish to store the Heat Pad elsewhere, store it in a dry place.
- Do not place heavy objects on top of the Heat Pad as it may cause the Heat Pad to crease badly, and damage the wiring.
- Do not use any moth-proofing chemicals on your Heat Pad.

After Sales Service

Your Elektra Comfort Heat Pad is not user-serviceable.

If the light does not illuminate and the Heat Pad does not heat up, check to see that:

1. You have followed the instructions properly.
2. The outlet supply is switched on.
3. The main supply fuse/circuit breaker has not blown/tripped.

If the appliance still does not work, contact the Service Centre.

If the light does illuminate and the Heat Pad does not heat up, check to see that:

1. The temperature control device is set correctly to setting 1, 2 or 3.

If the appliance still does not work, contact the Service Centre.

If the power supply cord or any wiring is damaged, or the Heat Pad has been misused or damaged, stop use immediately and contact the Service Centre.

| |
|---|
| TEL: 011 393 1328 FAX: 011 393 1501 EMAIL: technical@sp-africa.com |
|---|

Guarantee

Any defect affecting the functionality of the appliance which becomes apparent within one (1) year of its purchase, will be corrected by free repair or replacement, provided that the appliance has been used and maintained in accordance with the instructions, and has not been altered, abused or misused in anyway, and is accompanied with a Proof of Purchase to validate the age of the appliance.

Specifications

220 - 240V 50Hz 100W

The Environment

Please consider the environment when disposing of packaging.
Recycle wherever possible.

When you want to dispose of this appliance please do so responsibly, or contact us to assist with the disposal, or disposal advice.



**8 Palala Street, Protec Park, Terenure Ext. 42,
Kempton Park, South Africa, 1610**

Tel: 011 393 – 1328

Fax: 011 393 – 1501

www.sp-africa.com