

# ELEKTRA<sup>TM</sup>

---

## CARE



## INSTRUCTION MANUAL

PERSONAL SLIMLINE  
ELECTRONIC SCALE



# **Elektra Care Personal Slimline Electronic Scale**

Your new Elektra™ Care Electronic Scale has been designed to accurately measure your weight, and should provide many years of trouble free use if operated and maintained correctly.

Read these instructions before use and keep them safe. If you pass the appliance on, pass on the instructions too. Remove all packaging, but keep it until you are satisfied.

## **Caution**

- Do not use with wet feet.
- Do not use when the surface of the scale is wet.
- The surface of the scale is extremely slippery when wet.
- To avoid injury, ensure that the surface of the scale and your feet are dry before stepping on to the scale.

## **Specifications**

- Large LCD display
- Tempered glass
- Automatic zero re-setting
- Automatic switch-off
- Overload indicator
- Low battery indicator
- 1 x CR2032 Lithium battery, included.
- Maximum capacity – 150KG/330LB
- Division: 0.1kg/0.2lb

## **Power Supply**

- 1 x CR2032 Lithium battery, included.


Replace the batteries only when "**LO**" appears on the LCD display.

## **Directions for use**

- 1) Place the scale on a hard, flat surface (avoid carpeted floors).
- 2) To switch on the scale, tap the front centre of the platform with your foot and wait a few seconds until the zero appears on the display.
- 3) When the zero appears the scale is now ready for weighing.  
Step on the scale.
- 4) Stand still on the scale with your weight evenly distributed, and wait while the scale measures your weight.  
Your weight will then be displayed.

## IMPORTANT:

- 1) If you step on the scale **before** the zero appears on the display, the message "**Err**" will be displayed. If this happens, step off the scale and wait until it has switched itself off before repeating step two.

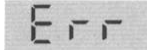

- 2) To ensure the most accurate reading:
  - (i) Always weigh yourself on the same scale placed on the same floor surface. Please note that uneven floors may affect readings.
  - (ii) Weigh yourself without clothing and footwear, before meals and always at the same time of day.

## Units of Measurement

Your Elektra Care scale can display your weight in **kilograms, pounds and stone**. The selector switch can be found on the underside of the scale. The smallest units of measurement are 0.1kg/0.2lb

## Automatic Switch Off

When you step off the scale it will automatically switch off after approximately 5 seconds.

## Overload Indicator

When the scale is overloaded, i.e. over its maximum capacity, the message "**O-Ld**" will appear on the screen.

## Low Battery Indicator

When the battery needs replacing, the message "**LO**" will appear on the screen.

To change the battery remove the battery cover on the back of the scale and carefully remove the battery from the battery compartment by drawing them sideways out from under the contacts.

Do not bend the contacts back whilst removing the battery as they may break. Please dispose of batteries in the correct manner.

## Care Instructions

- Clean your scale with a damp cloth but do not allow water to get inside the scale.
- The top of this scale becomes slippery when wet - please take care when weighing yourself.
- Do not use chemical cleaning agents.
- Do not allow your scale to become saturated with water as this can damage

the electronics.

- Keep your scale in a cool dry place.
- Keep your scale in the horizontal position when storing it.
- Remove the battery if the scale is not going to be used for a long period of time.
- Do not jump on or drop the scale.
- This product can only be used with static pressure. Avoid breaking the glass by striking, shocking or dropping.
- Children should only use the scale whilst under adult supervision.

## **After Sales Service**

The appliance is not user-serviceable. If it is not working properly check:

1. You have followed the instructions properly.
2. The battery is the correct battery and has charge.
3. The battery is making contact with terminal.

If the appliance still does not work, contact the Service Centre.

TEL: 011 393 1328 FAX: 011 393 1501 EMAIL: <a href="mailto:technical@sp-africa.com">technical@sp-africa.com</a>
---

## **Guarantee**

Any defect affecting the functionality of the appliance which becomes apparent within one (1) year of its purchase, will be corrected by free repair or replacement, provided that the appliance has been used and maintained in accordance with the instructions, and has not been altered/modified, abused or misused in anyway, and is accompanied with a Proof of Purchase to validate the age of the appliance.

## **Specification**

Model: 3100



**8 Palala Street, Protec Park, Terenure Ext. 42,  
Kempton Park, South Africa, 1610**

**Tel: 011 393 – 1328**

**Fax: 011 393 – 1501**

**Email: [info@sp-africa.com](mailto:info@sp-africa.com)**

**[www.sp-africa.com](http://www.sp-africa.com)**