

CARMEN

Classic



INSTRUCTION MANUAL

MASSAGE THERAPY FOOT SPA

Read these instructions before use and keep them safe. If you pass the appliance on, pass on the instructions too. Remove all packaging, but keep it until you are satisfied.

IMPORTANT SAFETY INSTRUCTIONS:

When using an electrical appliance, basic precautions should always be followed, including the following: Read all instruction before using this appliance.

DANGER -To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
2. Do not reach for an appliance that has fallen into water. Unplug it immediately.
3. Do not use while bathing or in a shower.
4. Do not place or store the appliance where it can fall or be pulled into a tub or sink. Do not place or drop into water or other liquid.

WARNING: -To reduce the risk of burns, fire, electric shock or injury to persons:

1. An appliance should NEVER be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
2. DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrical shock or injury to persons.
3. Close supervision is necessary when this appliance is used by, on, or near children, invalids or disabled persons.
4. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended; specifically any attachments not provided with the unit.
5. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the Service Center for examination and repair.
6. Do not carry this appliance by the supply cord or use the cord as a handle.
7. Keep the cord away from heated surfaces.
8. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair and other debris.
9. Never drop or insert any object into any opening.
10. To disconnect, turn all controls to the "0" (off) position, then remove the plug from outlet.
11. Do not use outdoors.
12. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
13. Do not attempt to plug or unplug unit while feet are in the water.
14. Do not stand on or in the foot massager. Use only while seated. Place feet in massager and remove feet from massager only while in a seated position.
15. Never operate on a soft surface, such as a bed or couch, where air openings may be blocked.
16. Connect the appliance to a properly grounded outlet only.
17. If supply cord is damaged, it must be replaced by the manufacturer or its service agent.
18. All servicing of this foot spa must be performed by authorized service personnel only.

19. Keep all packaging and plastic bags out of the reach of children as these pose the risk as a choking hazard.
20. The installation of a residual current operated protective device with a nominal trip current of 30mA offers additional protection. Contact your electrician.

FOR HOUSEHOLD USE ONLY

OPERATING INSTRUCTIONS

CAUTION

- The Massage Therapy Foot Spa should not be used if you have the following:

Tuberculosis, benign and malignant tumors, hemorrhages, inflammations of the skin, phlebitis and thrombosis, open and fresh wounds, bruises, broken skin or varicose veins, undetermined calf pain, or joint pain.

- Consult your physician if you have any therapeutic questions regarding the foot bath.
- Do not use if you are experiencing unexplained pain or swollen muscles, or following a serious muscular injury before consulting your physician.
- Keep away from children.
- Discontinue use of the foot bath if pain increases or swelling occurs after application.
- If you are pregnant or have a pacemaker, please consult your physician before using.


USER INSTRUCTIONS

1. Place foot bath on floor, fill foot bath with warm or cold water and plug into wall outlet. Never plug or unplug the foot bath while your feet are in water.
2. The MASSAGE THERAPY FOOT SPA can only be used with water.
3. Sit in a chair and place your feet in the water. Never stand in the foot bath.
4. Push the ON /OFF button once to switch the Foot Spa ON.
5. After use push the power button again to switch Foot Spa OFF and disconnect from outlet before emptying water out of the front of the unit.
6. Fill Footbath with warm water or cold water to the level indicator "Max". Do not fill above water-level indicator. Never overfill.
7. Pour water out of the footbath from the Easy-pour spout.

RECOMMENDED MASSAGING TIMES

Generally massage should not last longer than 20 minutes. Treatment may be applied once or twice daily as desired. Allow unit to cool before using again. As a person's metabolism varies, so does his or her response to massage. You will soon adjust your massage time to suit your needs and personal comfort. If you have any concerns regarding treatment, we recommend you consult your physician.

KNOW YOUR MASSAGE THERAPY FOOT SPA

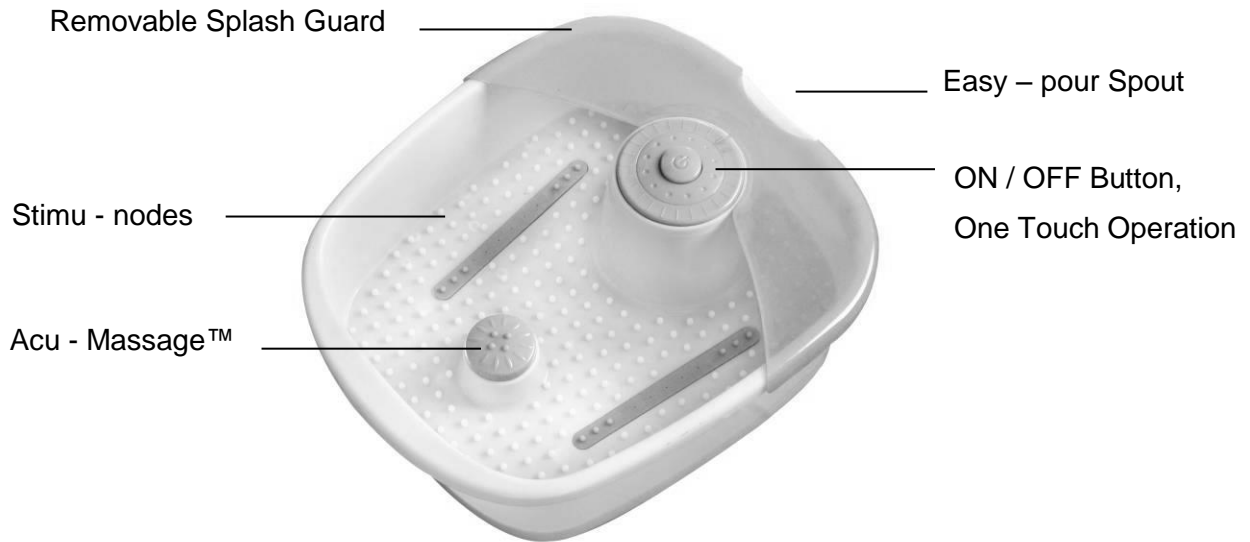

- **Stimu-nodes Massage**
 - Stimu-nodes to stimulate your feet increasing blood circulation.
- **Hydrotherapy - Bubbles**
 - Refresh tired feet.

Care and Maintenance

1. Switch the appliance off and unplug it.
2. Wipe the body of the appliance off with a soft damp cloth. If necessary use a little washing-up liquid.
3. Do not immerse the footbath in water or any other liquid.
4. Do not use harsh or abrasive cleaning agents or solvents.

After Sales Service

The appliance is not user-serviceable. If it is not working properly check:

1. You have followed the instructions properly.
2. The outlet supply is switched on.
3. The main supply fuse/circuit breaker has not blown/tripped.

If the appliance still does not work, contact the Service Centre.

Guarantee

Any defect affecting the functionality of the appliance which becomes apparent within one (1) year of its purchase, will be corrected by free repair or replacement, provided that the appliance has been used and maintained in accordance with the instructions, and has not been altered/modified, abused or misused in anyway, and is accompanied with a Proof of Purchase to validate the age of the appliance.

Specification

220 - 240V~ 50Hz 60W

The Environment

This appliance is rated at an operating input of 60watts.



Please consider the environment when disposing of packaging.

Recycle wherever possible.

When you want to dispose of this appliance please do so responsibly, or contact us to assist with the disposal, or disposal advice.



**8 Palala Street, Protec Park, Terenure Ext. 42,
Kempton Park, South Africa, 1610**

Tel: 011 393 – 1328

Fax: 011 393 – 1501

www.sp-africa.com